

C220
Coagulation Analyzer



**COMPACT DESIGN | ONE-CHANNEL
SCATTERED LIGHT | INTERNAL PRINTER**

- Easy keypad operation with LCD display Internal thermal-sensitive printer
- Reagent position with magnetic rotor
- Scattered light principle and percentage analysis ensures good results
- Open reagent with low consumption, close reagent on request
- Large memory to store up to 6000 test results
- QC and calibration program incorporated
- Auto power-on self-check
- Multi-language software available on request
- Electronically-linked pipette optional

C220

Coagulation Analyzer



Electronically-linked pipette




Internal printer



Keypad



Cuvette

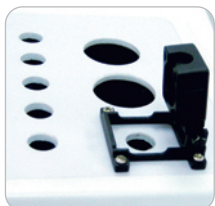

Channel-cover & Pipette-holder

Test Report			
Patient No.: 1	Name: Steven	Sex: Male	
Age: 32	Blood Type: AB	Bed No: 023	
No.	Item	Result	Range
01	PT (Prothrombin Time)	7.9s-0.64-12.45INR-21.9%	11.0 ~ 15.0
02	DFIB (Derived Fibrinogen)	7.9mg/dl	9.0 ~ 20.9
Test Date: 2008-06-05			

SPECIFICATIONS

Principle	Scattered light with percentage analysis
Test Channel	1 channel
Light Source	Durable LED detector
Test Time	Normally 20-60 secs, Maximum time up to 600 secs
Sample Position	5 positions
Reagent Position	2 positions (1 magnetic stirrer included)
Position Warming-up	37°C(±1.0°C)
Memory	6000 test results
Interface	RS-232
Display	LCD display
Input	Easy operation keypad
Output	Internal printer
Power Supply	AC 220V±22V,50Hz±1 Hz or AC 110V±22V,60Hz±1Hz
Net Weight	5KG
Dimensions L x W x H (mm)	270 x 270 x 120



Marketed By

Bioline Diagnostics LLP

H-1478, DSIDC, Narela Industrial Area, Delhi - 110 040

T. +91 11 4170 8100 | E. info@biolinediagnostics.com

www.biolinediagnostics.com