


FirstView®

One Step Rapid Test Device

Rapid test for HCV

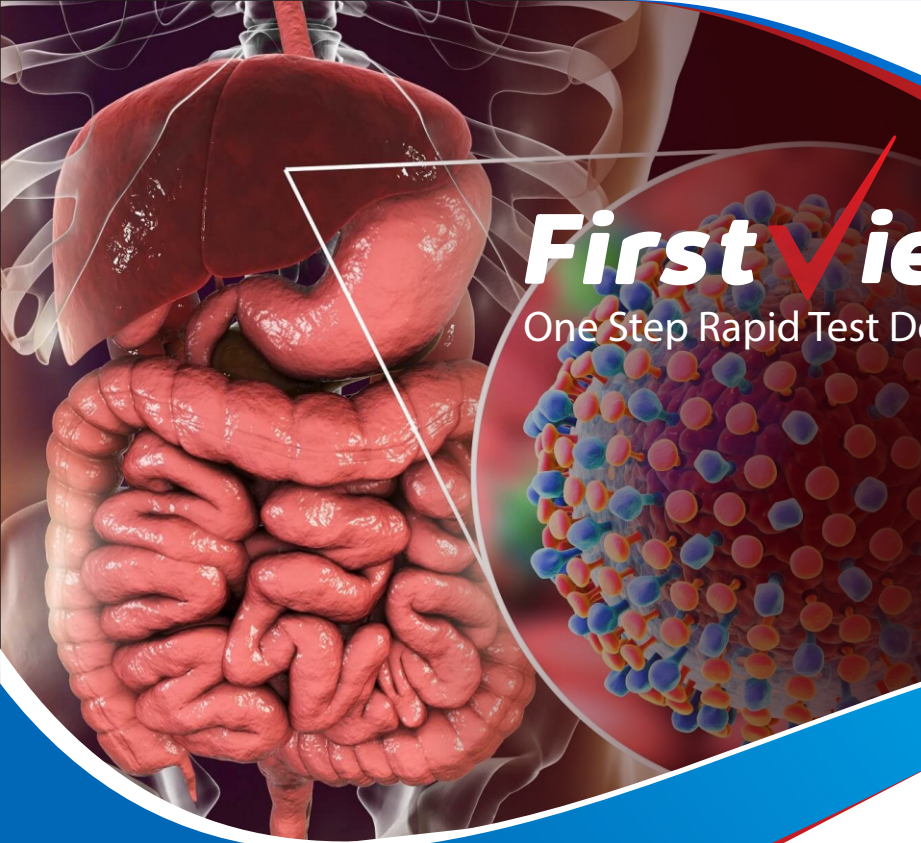


**For Qualitative
Detection of
HCV
Whole Blood
Rapid Test**

- Use of highly purified HCV antigens for Core, NS3, NS4, NS5 immobilized on the device.
- Run with PLASMA or SERUM
- Rapid test result in just 15 minutes
- Semi-quantitative point-of-care Test for HCV



READING TEST RESULTS

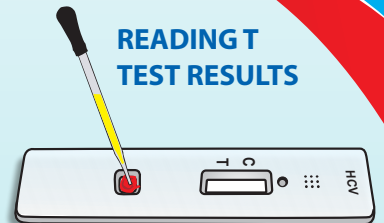


Add 1 drops (10 µL) of serum or plasma in sample well "S" using the dropper provided.



Add 2 drop of the Assay buffer



Result

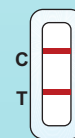

Positive



Negative



Invalid